

Ultrasonic Face Scrubber & Exfoliator



User Guide
NuLifeBeauty.com

Contents

Device Parts	4
How it Works	5
Charging	6
How to: Cleansing.....	7
How to: Mositurizing.....	8
How to: Lifting.....	9
Cleaning & Maintenance	10
Troubleshooting & Technical Issues.....	11
FAQ.....	12
Contact Us.....	13

Device Parts



Cleaning head allows for deep pore cleaning



One touch button controls varying modes



USB charging jack with covering



How it Works

There are 3 main functions to the Ultrasonic Face Scrubber & Exfoliator: Deep Skin Cleansing, Firming and Lifting, and Skin Rejuvenating. Each function is used differently and serves a different purpose.

- **Deep Skin Cleansing** - The high-frequency ultrasonic waves (Cleansing mode) will help open up your pores and clean them of any blackheads, dead skin, or dirt. Allowing your skin to regain its clearness and glow.
- **Firming and Lifting** - The EMS micro-current (Lifting mode) will gently stimulate facial muscles, tightening skin, lifting facial contours and reducing the appearance of fine lines and wrinkles.
- **Skin Rejuvenating** - The same ultrasonic waves (Moisturizing mode) that help open up the pores combined with negative ions help your skin absorb skincare products better and, promote collagen production and elastin in the skin, firming and tightening it up in the process.

Prior to using any of the three modes, you will need to ensure the device is charged. The device can be charged using the provided USB cable.

Charging

1 Plug charging cable into a USB port.

The battery will not be fully charged upon receiving the product. Charge the battery fully before use.

2 The USB can be plugged into a computer, adapter or charging block (sold separately).

LED light flickers continuously while charging. The light stops flashing when device is fully charged. Confirm the state of the battery before charging and tear the protective film down before first use.



- Device needs 30-60 minutes to fully charge. First time charging may take longer.
- Device can be used for up to 80 hours after full charge.
- Product will enter standby mode after 10 mins and turn off after 1 minute if not used.

How to: Cleansing

- 1** Press and hold power for 2 seconds to turn on device.
- 2** For cleansing mode click the power button 1 more time.
- 3** Ensure face is clear of makeup, apply cleanser to face or ensure face is wet.
- 4** Scrub face in upward motion with device head.



How to: Moisturizing

- 1** Press and hold power for 2 seconds to turn on device.
- 2** For moisturizing mode click the power button 2 more times.
- 3** Ensure face is clear of makeup, then apply thin layer of skincare product for absorption.
- 4** Press Sonic Head of device around face.



How to: Lifting

- 1** Press and hold power for 2 seconds to turn on device.
- 2** For lifting mode click the power button 3 more times.
- 3** Ensure face is clear of makeup, then gently press Sonic Head in an upward lifting motion.




Cleaning & Maintenance

- 1** Make sure device is powered off.
- 2** Clean device with a cotton pad or soft fabric.
Do not use solvent or any corrosive product to clean Sonic Head to avoid damage.
- 3** Keep device stored in cool, dry place.
Do not store the device in a wet place such as in the bathroom. Keep the device in a location away from children.



Troubleshooting and Technical Issues

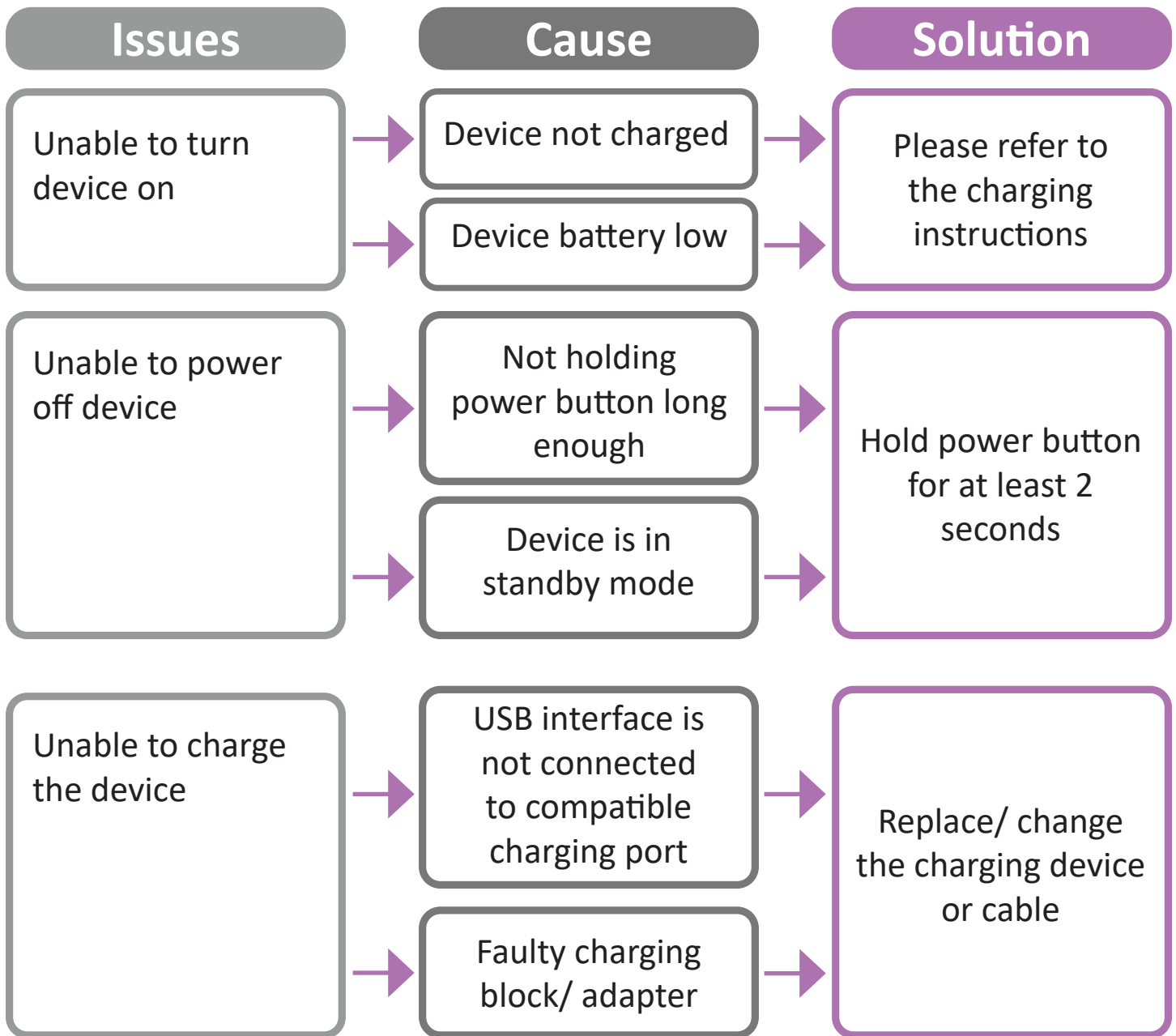

FAQ

Do I need to use a face product with the device?

The device is specifically used to enhance moisturizing and cleansing routines. It is safe to use with a thin layer of cleanser, moisturizer or water.

How do I know if it is working?

You will notice the accumulation of pore residue on the Sonic Head after use. The UV light below the Sonic Head highlights the residue.

How often can I use the device?

It is recommended to use the device once a day for best results.

Is it supposed to get hot?

The device may get warm as you are using it. It shouldn't be hot or uncomfortable on the skin.

Are there any side effects?

There are virtually no side effects. In rare cases, some users may experience skin redness and irritation. If that occurs please discontinue use or reduce frequency.

Contact Us

If you have any questions or concerns, please do not hesitate to reach out to us!

We will be happy to help!

Email us:

help@nulifebeauty.com

Message us on Facebook or Instagram:

@NuLifeBeautyUSA

